

FIRST-CALL MEDICAL, INC.

The call that saves...


574 Boston Rd #11
Billerica MA 01821
Phone: 800-274-5399
Fax: 978-670-5457
Website: www.EKGreport.com

HEARTRAK SMART Patient Manual

TO RECORD YOUR ECG

- Be sure that the skin on your chest is free of any oils, lotions, or talc powder before applying electrode patches.
- The WHITE tip lead wire is placed beneath the right collarbone as in the illustration to the right.
- The BLACK tip lead wire is placed approximately 2 inches below the left breast area as in the illustration to the right.
- You can record ECG event at any time by pressing and releasing the white RECORD button found on the top of the recorder. The recorder will produce a tone with the rhythm sound of your heart beat.
- Please remain as still as possible until the tone stops. The recorder captures ECG before the RECORD button is pressed and after the RECORD button is pressed.
- After you RECORD your ECG, the recorder automatically returns to scanning mode. You should try to call the RECEIVING CENTER after each recorded event.


Electrodes and Lead Wire Placement

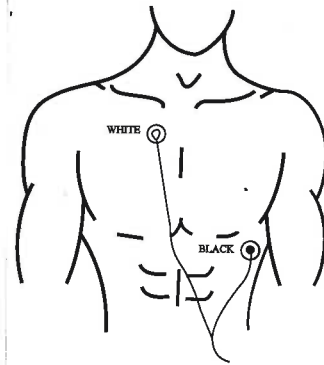

TO TRANSMIT YOUR EKG

- Call the Receiving Center at 1-800-274-5399. The Center is open 24-hours a day, 7 days a week.
- Please identify yourself to the technician and explain the symptom(s) that you recorded.
- Remove the monitor from the carrying case.
- Upon instruction from the technician, place the telephone mouthpiece (the part you talk into) over the small speaker hole found on the front of the recorder.
- Press and release the SEND button. The recorder will produce a steady tone. It is now playing back your ECG.
- When the tone stops, please pickup your telephone handset to speak with the technician.
- If the ECG signal was received successfully, you will be instructed to press the RESET button to clear the recorder's memory.

Mouthpiece of Phone on Recorder

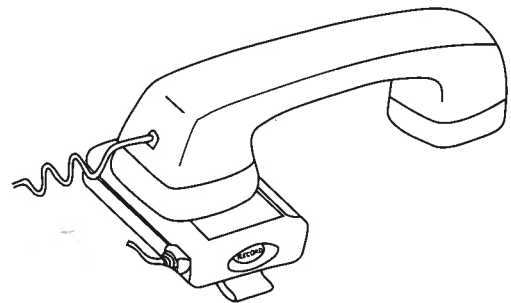

BATTERY REPLACEMENT

If the BATTERY LOW light is lit on the front of the recorder, please replace the battery by sliding the panel on the side of the recorder. Insert a new "AA" battery according to the polarity guide inside.

Sliding the Battery Door Panel

