

# FIRST-CALL MEDICAL, INC.

The call that saves...


**574 Boston Rd #11**  
**Billerica MA 01821**  
**Phone: 800-274-5399**  
**Fax: 978-670-5457**  
**Website: [www.EKGreport.com](http://www.EKGreport.com)**

## HEARTRAK SMART Patient Manual

### TO RECORD YOUR ECG

- Be sure that the skin on your chest is free of any oils, lotions, or talc powder before applying electrode patches.
- The WHITE tip lead wire is placed beneath the right collarbone as in the illustration to the right.
- The BLACK tip lead wire is placed approximately 2 inches below the left breast area as in the illustration to the right.
- You can record ECG event at any time by pressing and releasing the white RECORD button found on the top of the recorder. The recorder will produce a tone with the rhythm sound of your heart beat.
- Please remain as still as possible until the tone stops. The recorder captures ECG before the RECORD button is pressed and after the RECORD button is pressed.
- After you RECORD your ECG, the recorder automatically returns to scanning mode. You should try to call the RECEIVING CENTER after each recorded event.


### Electrodes and Lead Wire Placement

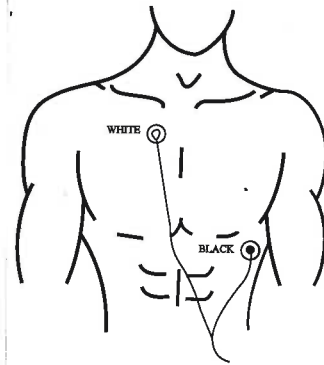

### TO TRANSMIT YOUR EKG

- Call the Receiving Center at 1-800-274-5399. The Center is open 24-hours a day, 7 days a week.
- Please identify yourself to the technician and explain the symptom(s) that you recorded.
- Remove the monitor from the carrying case.
- Upon instruction from the technician, place the telephone mouthpiece (the part you talk into) over the small speaker hole found on the front of the recorder.
- Press and release the SEND button. The recorder will produce a steady tone. It is now playing back your ECG.
- When the tone stops, please pickup your telephone handset to speak with the technician.
- If the ECG signal was received successfully, you will be instructed to press the RESET button to clear the recorder's memory.

### Mouthpiece of Phone on Recorder

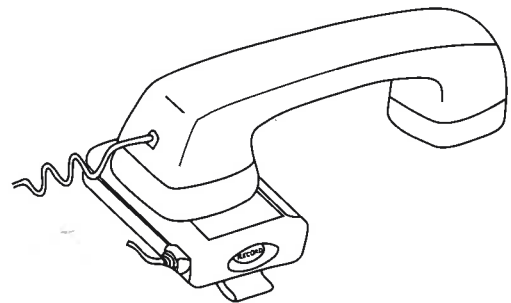

### BATTERY REPLACEMENT

If the BATTERY LOW light is lit on the front of the recorder, please replace the battery by sliding the panel on the side of the recorder. Insert a new "AA" battery according to the polarity guide inside.

### Sliding the Battery Door Panel

